


Expandable Bariatric

Extra high SWL expandable alternating pressure mattress for the prevention and treatment of pressure injuries.

- Expandable width 36" / 42" / 48"
- Expandable length 80" / 85" / 90"

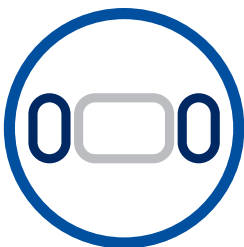

SIDE BOLSTERS

Increased lateral support and stability with reinforced side bolsters



WEIGHT RANGE

Superior weight range capacity up to 454 kg



AUTOFWOWER

Automatically boost pressure to protect the sacral area during seated position



QUICK RELEASE CPR

Easy opening and accessible CPR unit for rapid deflation

SAFETY CARE

Foam underlay providing safety from bottoming out

EXPANDABLE

Expandable length and width for transport

NYLON CELLS

Designed for durability and increased weight capacity

SPECIFICATIONS

CODE	APMPC-B01S
SIZE	Expandable Bariatric
CAPACITY (SAFE WORKING LOAD)	40-454 kg
SYSTEM ALTERNATION	1 in 2
WARRANTY	1 year
PRESSURE INJURY RISK	High to very high risk
CERTIFICATION	IEC 60601, CE Mark
ARTG	370834
MATTRESS DIMENSIONS (L x W x H)	2032-2286 x 914-1219 x 250 mm

