


AreaCare Comfort

A softer three layer foam surface for clients with specific comfort needs. Featuring a more immersive top layer memory foam with an updated support layer available in Single and King Single models.

	AREACARE COMFORT SINGLE	AREACARE COMFORT KING SINGLE
LENGTH	1960 mm	1960 mm
WIDTH	865 mm	1050 mm
HEIGHT	150 mm	150 mm
LAYERS	3	3
SIDE WALLS	✓	✓
SWL	200 kg	220 kg
INFECTION CONTROL SEAMS	✓	✓
HINGED		✓
CODE	BSPS4003	BSPS4002
WARRANTY	6 year body; 2 year cover	

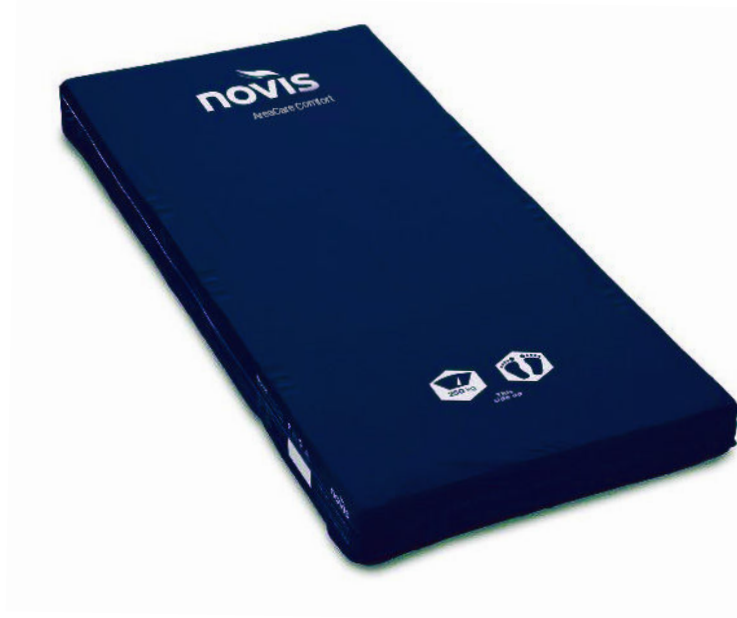

Top Cover

Latex free, multi-stretch polyurethane coated top cover for reduced shear and friction.

Top Layer

Increased comfort memory foam top layer allows immersion

Support Layer

High resilience foam for maximum pressure redistribution.

Core Layer

Firm core layer and edges provide form, strength and stability

Hinged Base

Six hinge system to conform to profiling beds. (King Single only)

