

Prodigy® Mattress Overlay



Drainage holes in
base to draw away fluids

Non-powered for easy portability

Durable neoprene construction for easy cleaning

The ROHO Prodigy Mattress Overlay is an adjustable, zoned, air filled pressure reducing surface suitable for patients at low to medium risk.

- 3 independently adjustable, interconnected zones
- durable polyvinyl construction for easy cleaning with soap and water
- non-powered for easy portability
- low friction and shear environment
- lightweight overlay (0.9 kg per section; 2.7 kg complete system)
- drainage holes in base to draw away fluids
- suitable for users to maximum capacity of 135 kg



ROHO's air floatation technology redistributes forces that can cause peak pressure, shear and friction, providing effective skin protection and relief for patients who are at risk of developing pressure or tissue injuries.

The Prodigy Mattress Overlay system is comprised of three interconnected sections (each 680 x 890 mm) made from independent, polyvinyl air cells that encourage air movement and ventilation. The sections easily snap together to form an overlay. The 90 mm high cells can be adjusted per segment to address the specific needs of the individual patient body size and shape.


The durable polyvinyl construction is flame-resistant and easy to clean with soap and water. Drainage holes in the overlay base help to pull fluids away from incontinent patients.

Prodigy® Mattress Overlay

Optional
Enclosure Cover



3 independently adjustable,
interconnected zones

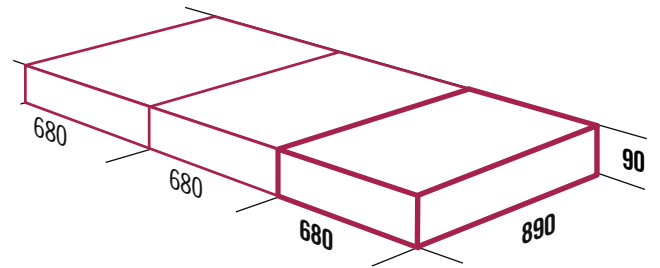

Includes a hand inflation pump, repair kit and operating instructions. The ROHO Mattress Overlay Enclosure Cover is also available as an optional accessory – a fluid resistant, fully enclosed mattress overlay cover for increased infection control and incontinence protection. Straps secure the cover to a standard single bed mattress, while carry handles help with product transport.



ITEM	Mattress Overlay	Enclosure Cover
CODE	R-PRODIGY	R-COV-MR2PCSTD-P
CAPACITY	135 kg	
LENGTH	2075 mm	
WIDTH	915 mm	
HEIGHT	90 mm	
WEIGHT	4.3 kg	
CONFORMS TO	ARTG 111257	
WARRANTY	1 year	

